



## INSTALLING CEREBRUM NEXUS 2.1.3

### NOTICE

Cerebrum assumes no liability whatsoever, and disclaims any warranties whether express or implied, written, oral, statutory or otherwise relating to the information, recommendations or technical feedback provided herein and its use, including any liability for warranties relating to fitness for a particular purpose, performance, quality, merchantability, or non-infringement of any intellectual property right. The user is responsible for determining the suitability of Cerebrum information, recommendations or technical feedback for the user's products and intended applications.

## Table Of Contents

Table Of Contents.....	2
Document Information.....	3
Change History .....	3
About This Guide.....	4
Contacting Support .....	4
Installation Guide.....	4

## Document Information

### Change History

The following table summarizes major changes to this document.

Document Version	Date	Major Changes
1.0	2022/06/29	Initial version
1.1	2022/08/29	Updating to Nexus 2.0.1
1.2	2023/01/26	Updating to Nexus 2.1.0
1.3	2023/02/13	Updating to Nexus 2.1.1
1.4	2023/02/16	Updating to Nexus 2.1.2
1.5	2023/02/23	Updating to Nexus 2.1.3

## About This Guide

This guide was created to assist you install Cerebrum Nexus 2.1.3. This version of our software introduces many new features and primarily introduces a new security handshake between Cerebrum Nexus and LABdivus, and it looks for and installs Oracle Java if needed. We recommend that you use this guide to plan and execute the upgrade and contact Cerebrum Support with questions before or during the upgrade.

## Contacting Support

Online:

<https://support.cerebrumcorp.com>


Call:

(855) 469-5467

Mon-Fri 7am-7pm

## Installation Guide

1. Please make sure you are an administrator on the computer or can get administrative rights. You'll also need to be an admin on LABdivus (Lab Manager or Pathologist).
2. Click on Start and type in "add or remove" and the Add or Remove Programs Settings should appear. Click on it. Scroll down to Cerebrum Nexus and Uninstall

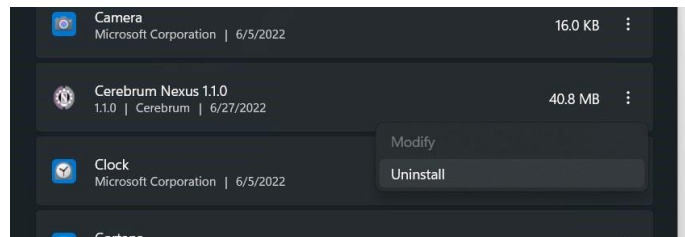

3. Download Cerebrum Nexus 2.1.3 from <https://www.cerebrumcorp.com/download-nexus>.


[Download Cerebrum Nexus 2.1.3](#)

[Download the 2.1.3 Install Guide](#)

Released 02-23-2023


4. Double click on the downloaded file to start the installation. You'll be warned by Windows Smart Defender that the program is not recognized by Windows. Click "more info" and then run the program anyway.

**NOTE:** If you don't have Oracle Java installed, Nexus will download it automatically and start the installation, you'll need to click on those prompts from that installer. Then Nexus will automatically resume installation.

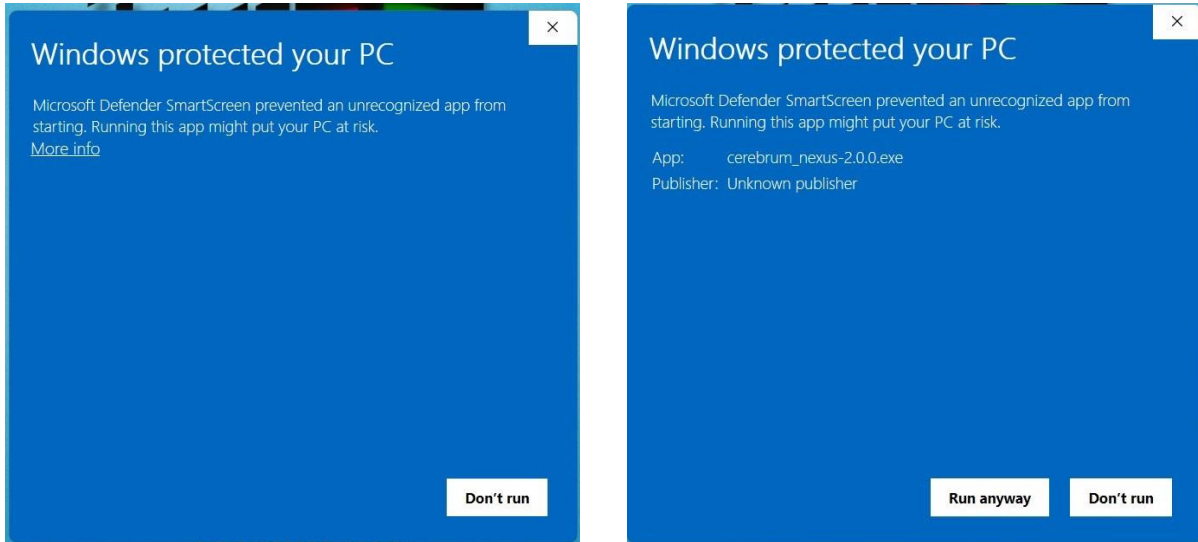

5. After the installation, you'll need to run Nexus by clicking on "Start" and searching for nexus, then clicking on the icon.
6. You'll be prompted with a pop up asking for your LABdivus Key. You can find this buy logging in and navigating to Configuration > Application > Global Variables

Ldap Host	<input type="text"/>	LDAP Host
Ldap Port	<input type="text" value="389"/>	LDAP Port
LABdivus key	<input type="text" value="aMz0FfuZM+o95eMXfni84dsxBRZD7Y05OKnVx15zwwISOOnLUfxLh8Cr"/>	Random Nexus Instance Key
Package Format	<input type="text" value="zip"/>	Package Format (tar/zip)
Package Limit	<input type="text" value="20"/>	Max Amount of Files in Package

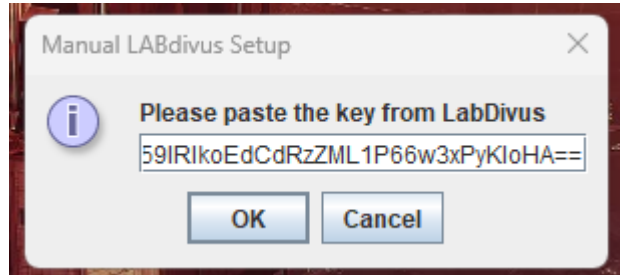

7. Close your web browser completely, reopen it, and LABdivus to complete the process.